

Hand-held beauty apparatus



装
订
线

Table of contents

1. Product overview and parameters.....	01
2. System interface.....	02
3. Operation steps.....	02
4.The operation method.....	05
5.The use of taboos.....	07
6.Matters needing attention.....	08
7. Instrument repair and maintenance.....	09
8. Warranty card.....	10

1. Product overview and parameters

Welcome to the hand-held beauty apparatus, the hand-held beauty apparatus is a dual-use beauty instrument for home and home. Please read this manual carefully before use and operate this beauty device strictly according to the instructions in the manual.




Product name: Hand-held beauty apparatus	Power: 28-56W
Product size: 21*42*46cm	Product model: z-1808
Power supply voltage: 220-240V	Working voltage: 28V
Operating head: 2.0/1; 3.0/1; 4.5/1*15,000 rounds per operating head	
Implementation standard:	
GB/T 36419-2018 (Household and Similar Skin Beauty Devices)	
GB 4706.1-2005 (Safety of Household and Similar Electrical Appliances Part 1: General Requirements)	
GB 4706.15-2008 (Special requirements for safe skin and hair care appliances for household and similar electrical appliances)	

2. System interface



3. Operation steps

-  Clean the operating area first (if you have makeup, use makeup remover to remove makeup first, then use facial cleanser to clean it);
- There are three operating heads, 2.0/3.0/4.5 respectively. Connect the operating heads to be used to the instrument according to the hole positions. At the beginning of operating the face, first choose 4.5mm depth to operate;

3. Connect the power cord, press the circular button above the power cord to turn it on (two circular buttons on the instrument, the one with up and down arrows is the gear button (⏮), and the one with only one triangle pattern is the start button (⏻)), press the gear to adjust Press the button to turn on the instrument, if you need to shut down, please press and hold the gear button for 2 seconds (⏮);

4. The amount of energy can be adjusted after starting up (⏮), there are 9 gears in total. It is recommended to use the lowest gear for testing operations at the beginning; after pressing the start button of the instrument, the instrument is working and the energy is outputting, and the instrument will beep once a second.;

5. Test the energy intensity: first pull the right lower jaw at a constant speed to the top of the corner of the mouth. If there is no feeling to increase the energy, it is normal to have a feeling of heat, tingling or soreness during operation. Lower the energy (adjust the amount of energy according to personal conditions), because the position with more bones is more sensitive, and the energy received will be relatively low (different operating parts, the amount of energy received is also different, you need to adjust the amount of energy according to your own situation.);

Normal skin: For the first use, it is recommended to use the energy at a low level first, and adjust the energy to the maximum that the individual can bear according to the feeling during the operation;

Red blood skin, sensitive skin with thin stratum corneum: the first time to use the energy to operate at a low level; after 2-3 operations, adjust the energy to the maximum that the individual can bear according to the situation of the individual skin;

4. The operation method





6. During the operation, it is necessary to apply the medium on the operation site. During the operation, the probe must be close to the skin, and the operation should be carried out by moving up and down at a uniform speed. After the medium is absorbed or dried during the operation, it needs to be replenished immediately;

7. The normal operation time of the full face is about 30-45 minutes; (the size of each person's face is different, so the operation time is also different);

8. After using the instrument, you need to use a wet towel to clean up the residual medium on the probe;

9. After the operation, clean the face with water or lukewarm water, immediately put on a moisturizing mask, leave it on for 10-15 minutes, and then remove it. If the water on the mask paper is absorbed and dried in the middle, replenish the water in time; (not It is recommended to wash your face with hot water)

Note: If you need to operate other parts, please consult a professional before operation.

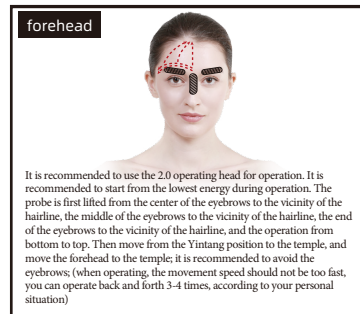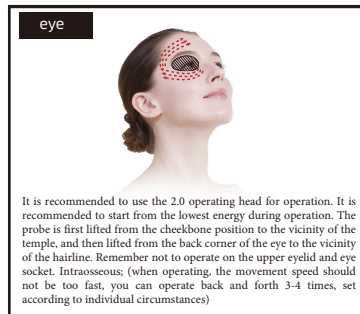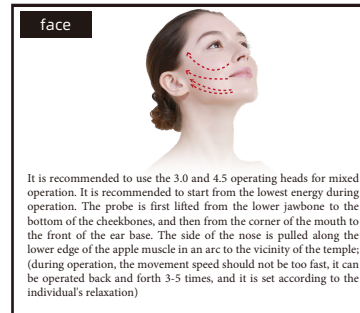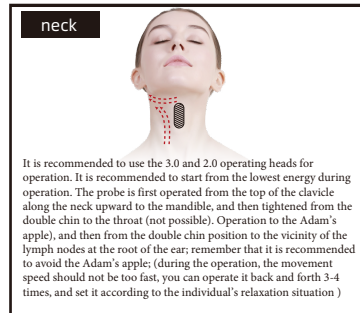

5.The use of taboos

▪ **Note:** During operation, the probe cannot be positioned and repeated at one position, and it is operated in a uniform movement mode. It is not recommended to adjust the energy to the highest during initial operation, it is recommended to start from the lowest energy. Then adjust the amount of energy according to the situation of the operating part, and adjust it to the energy range that the individual can accept; you will feel a slight heat sensation during the operation, and there will be a small tingling feeling in the more sensitive parts, which may be felt after the operation. A feeling of soreness or numbness is normal;

▪ For normal skin, it is recommended to use it once every 5 days for the first 5 operations. After 5 operations, it is recommended to use it once every 10 days. 10 times is a cycle. For sensitive skin, it is recommended to operate once every 10 days for the first 5 operations. After 5 operations, it is recommended to operate once every 15 days. 10 operations are a cycle. The number of operation cycles is set according to the degree of individual skin sagging and sagging;

▪ After the operation, different degrees of detoxification and metabolism may occur on the face: acne, pimples, redness, itching, slight edema, soreness, etc. are all normal, and normal skin generally repairs automatically within 1-3 days. Individual special cases will be repaired for a long time. In the above cases, you can apply more moisturizing face to the face. These are usually normal, mild and temporary reactions.

People with the following symptoms should use this instrument with caution. Please consult a doctor or professional before use. The details are as follows:

1. Clean the operating part before use to avoid dirt or dirt. Before operating, please remove the metal objects worn on the head to avoid unforeseen situations and affect the effect;
2. Menstrual period is not recommended;
3. It is not recommended to use when the face is allergic, it is recommended to use it after a week after the allergy is cured;
4. It is recommended to use the beauty device after 20 days of facial injection of hyaluronic acid or face-lifting needle;
5. It is recommended to avoid this area where the prosthesis and silicone are applied to the operation site;
6. It is not recommended to use the area with wounds or large pimples or pustules on the face;
7. Patients with serious diseases, heart disease, diabetes, mental illness, etc. should be used with caution
8. Patients with severe skin disease or transfection disease or hormone face should be used with caution;
9. It is not recommended to use it during pregnancy, and it can be used after the body of the confinement child recovers after childbirth;
10. It is not recommended to use gold wire (non-absorbable metal wire) implantation on the face;

6. Matters needing attention

- It is recommended to apply a hydrating mask at least 3-4 times a week, and drink no less than 2000cc (ie, eight glasses of water) every day;
 - Eat less spicy food, less smoking, drinking, staying up late, eat more fruits and vegetables, and eat lighter food;
 - It is not recommended to eat a large amount of heavy-pigmented things such as soy sauce within a week;
 - Note that it is not recommended to go to a high temperature environment such as a high temperature sauna and do not expose to the sun within 1 week after the operation;
 - Do not wash your face with hot water within three days (water not exceeding body temperature is sufficient);
 - No special care is required, just hydrate and sunscreen normally;
 - During the operation, the operation site should be kept moist, and the probe should be kept close to the skin;
 - When using on the face, it is not recommended to operate on the orbital bone, upper eyelid, temple or Adam's apple, bridge of the nose, etc.;
 - When operating, it is recommended to slowly adjust the energy from small to large, adjust to the most suitable energy intensity to operate, do not directly turn on the maximum energy operation, directly turn on the maximum energy operation, it may damage the skin (scalds).);
 - The use of this product may be harmful to people with allergic constitution, people with skin diseases and allergic dermatitis, and sensitive skin. The crowd, etc. have discomfort or adverse effects. If there are adverse effects during use, the appliance should be stopped immediately;
 - Do not use the device on areas with compromised skin integrity (eg burns, blisters, cuts, wounds, active skin disease, and/or infection);
 - This appliance should not be used on children.

WARNING: Prolonged or repeated use of this device on a surface may result in potential damage to the skin;

WARNING: Prolonged or repeated use of this device on a surface may result in potential damage to the skin;

WARNING: Do not bypass or cancel the appliance's inherent safety features - contains information on cleaning and disinfection (if necessary) of the appliance or its parts.

7. Instrument repair and maintenance

1. After each use, please clean the instrument with clean water and keep it properly.
2. Before using the instrument, a plug with a grounding pin must be used, and make sure that the power socket of the instrument is truly grounded.
3. Make sure that the voltage of the instrument is suitable. If the local power supply voltage is unstable, we recommend that the user add a regulated power supply with matching power.
4. In order to ensure the effect of the instrument and ensure the normal service life, please use the designated accessories provided or recommended by the original manufacturer.
5. The instrument should not be placed in a humid place or near water sources, and do not expose the instrument directly to sunlight.
6. Do not place the instrument close to a strong heat source, as this may affect the life and normal use of the instrument.
7. Media such as water, oil, etc. shall not penetrate into the interior of the instrument, and shall not be knocked, knocked or dropped heavily to avoid damage to the instrument.
8. The warranty period of the instrument manufacturer is one year. During the warranty period, the machine has non-artificial damage, and the manufacturer is responsible for free maintenance and replacement of parts (the return shipping needs to be borne by the manufacturer, and the replacement parts and maintenance costs are paid by the factory).
9. After the warranty period ends, the manufacturer will provide lifetime maintenance for the instrument (the manufacturer for replacement parts only charges the cost, not the maintenance fee).
10. It is forbidden to use things containing a lot of alcohol for wiping cleaning.

8. Warranty card

Warranty Card		
Product name:		model:
Purchasing date:		Store name:
Address:		
Customer name:	Tel:	wayyanty period:
Failure cause:		
Repair content:		
Warranty terms:		
1. If the product has quality problems caused by non-human damage, it will not be replaced if there is no quality problem.		
2. If it is a software failure, the machine can be repaired within three days without replacement.		
3. If the product has hardware failure caused by non-human factors, the company will provide warranty within one year from the date of purchase of the product.		
4. The machine does not turn on, no warranty.		
Remark:		